

Primerdesign™ Ltd

# Sandfly fever Sicilian virus

genesig® Easy Kit  
for use on the genesig® q16

50 reaction





GENESIG

Kits by Primerdesign





For general laboratory and research use only

# genesig® Easy: at a glance guide





## For each RNA test

Component	Volume	Lab-in-a-box pipette	
SFSV primer/probe mix	5 µl		
Your RNA sample	5 µl		
oasig OneStep Master Mix	10 µl		

## For each positive control

Component	Volume	Lab-in-a-box pipette	
SFSV primer/probe mix	5 µl		
<u>Positive control template</u>	5 µl		
oasig OneStep Master Mix	10 µl		

## For each negative control

Component	Volume	Lab-in-a-box pipette	
SFSV primer/probe mix	5 µl		
<u>Water</u>	5 µl		
oasig OneStep Master Mix	10 µl		

# Kit Contents



- **SFSV specific primer/probe mix (BROWN)**  
Once resuspended the kits should remain at -20°C until ready to use.



- **Lyophilised oasisig™ OneStep Master Mix**



- **Lyophilised oasisig™ OneStep Master Mix resuspension buffer (BLUE lid)**



- **SFSV positive control template (RED lid)**



- **Internal extraction control RNA (BLUE lid)**



- **RNase/DNase free water (WHITE lid)**



- **Template preparation buffer (YELLOW lid)**



- **54 x genesig® q16 reaction tubes**

## Reagents and equipment to be supplied by the user

### genesig® q16 instrument

### genesig® Easy Extraction Kit


This kit is designed to work well with all processes that yield high quality RNA and DNA but the genesig Easy extraction method is recommended for ease of use.

### genesig® Lab-In-A-Box

The genesig Lab-In-A-Box contains all of the pipettes, tips and racks that you will need to use a genesig Easy kit. Alternatively if you already have these components and equipment these can be used instead.


# Step-by-step guide

## 1. Resuspend the test components

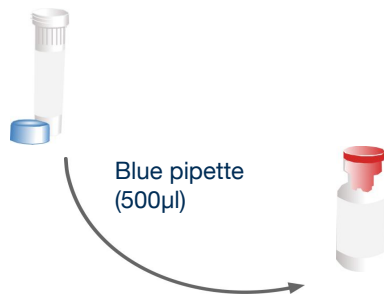

Use the blue pipette to transfer 500µl\* of the oasig OneStep Master Mix resuspension buffer into the tube of lyophilised oasig OneStep Master Mix and mix well by gently swirling.

\*Transferring 525µl of the oasig OneStep Master Mix resuspension buffer to your oasig OneStep Master Mix (instead of the 500µl recommended above) will enable you to take full advantage of the 50 reactions by accounting for volume losses during pipetting. In order to do so with the genesisig Easy fixed volume pipettes use 1x blue, 2x red and 1x grey pipettes to make the total volume. Please be assured that this will not adversely affect the efficiency of the test.

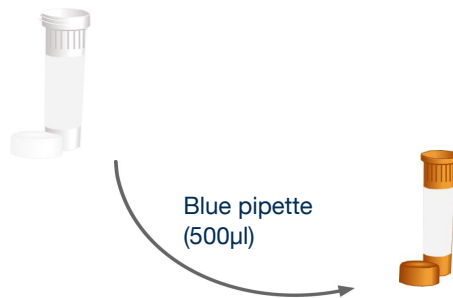

Then use the blue pipette to transfer 500µl of water into the brown tube labelled SFSV primers/probe. Cap and shake tube to mix. A thorough shake is essential to ensure that all components are resuspended. **Failure to mix well can produce poor kit performance.**

These components are now ready to use.

Store them in the freezer from hereon.

### Top tip

- Ensure that the primer/probe mix is mixed thoroughly before each use by shaking.
- Once resuspended do not expose genesisig Easy kit to temperatures above -20°C for longer than 30 minutes at a time.

## 2. Internal extraction control



Use the blue pipette to transfer 1000µl (2 x 500µl) of template preparation buffer into the Internal Extraction Control RNA tube. Cap and shake tube to mix.

Your kit contains Internal Extraction Control RNA. This is added to your biological sample at the beginning of the RNA extraction process. It is extracted along with the RNA from your target of interest. The q16 will detect the presence of this Internal Extraction Control RNA at the same time as your target. This is the ideal way to show that your RNA extraction process has been successful.

### **If you are using an alternative extraction kit:**

Use the red pipette to transfer 10µl of Internal Extraction Control RNA to your sample **after** the lysis buffer has been added then follow the rest of the extraction protocol.

### **If using samples that have already been extracted:**

Use the grey pipette to transfer 5µl of Internal Extraction Control RNA to your extracted sample.

## 3. Add primer/probe mix to all reaction tubes



For every reaction to be run, use the grey pipette to add 5µl of your SFSV primers/probe mix to every tube.

### Top tip

- Always pipette the primer/probe mix directly into the bottom of the tube.
- You can label the tube lids to aid your reaction setup but avoid labelling tube sides.

## 4. Add Master Mix to all reaction tubes



For every reaction to be run, use the red pipette to add 10µl of the oasis OneStep Master Mix to the tubes containing primer/probe mix.

Move swiftly to begin your q16 run, as any delay after the oasis OneStep Master Mix has been added can effect the sensitivity of your test.

### Top tip

- Always add the oasis OneStep Master Mix to the side of the tube to reduce the introduction of bubbles.

## 5. Negative control



For each test you will require a negative control. Instead of RNA water is used. This sample should typically prove negative thus proving that all of your positive samples really are positive.

To create a negative control reaction simply use the grey pipette to add 5µl of the water to the required reaction tubes. Close these tubes after adding the water.

Because some genesig kit targets are common in the environment you may occasionally see a “late” signal in the negative control. The q16 software will take this into account accordingly.

### Top tip

- Always add the water to the side of the tube to reduce the introduction of bubbles.

## 6. Set up a test




For each sample you wish to analyse, use the grey pipette to add 5µl of your RNA sample to the required reaction tubes. Close these tubes after adding the sample. Always change pipette tips between samples.

### Top tip

- Always add the RNA sample to the side of the tube to reduce the introduction of bubbles.

## 7. Positive control

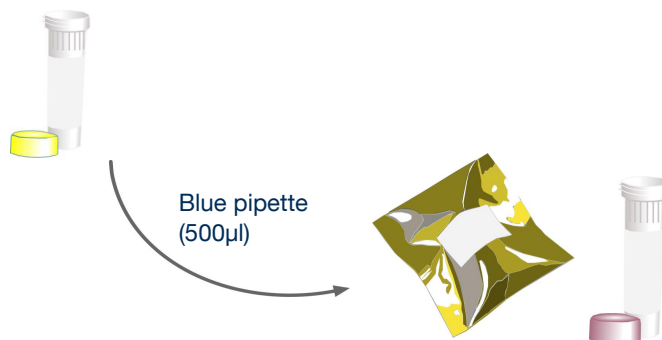

Use the blue pipette to transfer 500µl of template preparation buffer into the positive control template tube. Cap and shake tube to mix.

Each time you run a test you will require a positive control. This is a small portion of RNA from your target of interest. It serves two purposes:

1. It will always test positive so it shows that everything is working as it should be.
2. The q16 software knows how much RNA is present in the positive control. So it can automatically compare your sample of interest with the positive control to calculate the amount of target RNA in your sample.

To create a positive control reaction simply use 5µl of the positive control instead of your RNA sample.



Take great care when setting up your positive control. The positive control template has the potential to give you a false positive signal in your other samples. Set positive controls up last after all other sample tubes are closed. Always change pipette tips between samples. You may even choose to set up positive controls in a separate room.

### Top tip

- Always add the positive control template to the side of the tube to reduce the introduction of bubbles.

## 8. Running the test

Place the tubes into the correct positions in your q16 as defined by the software, this may include positioning of empty tubes to ensure that the q16 lid is balanced. The run can then be started.



### Top tip

- Before loading tubes into the q16, check for bubbles! Flick the bottom of the tubes to remove any bubbles that may have formed during the test setup.
- Apply centrifugal force with a sharp wrist action to ensure all solution is at the bottom of the reaction tube.
- When repeating a test you can use a previous file as a template by clicking 'open' then selecting File name > Files of Type > Experiment file as template



# What do my results mean?

Analysis of your data is carried out automatically by the genesig q16. The following information is designed to help you fully understand a result or to troubleshoot:

## “Positive”

### **Explanation**

Your sample has produced a positive result. Your target of interest is present and you can use the reported quantity.

## “Negative”

### **Explanation**

Your sample has produced a negative result. The target is not present in your sample.

## “Test contaminated”

### **Explanation**

The Negative Control should be completely free of any DNA/RNA. If you see this error message it means that at some point during the setup, the Negative Control has been contaminated with DNA/RNA and has given a positive signal. This contamination has invalidated the test. The Positive Control and your test samples are both possible sources of contaminating DNA/RNA. The genesig q16 reaction tubes from previous runs will also contain very high amounts of DNA so it is important that these are carefully disposed of after the run is completed and NEVER OPENED. It may be the case that your kits have become contaminated which will lead to the same problem occurring repeatedly.

### **Solutions**

1. Clean your working area using a commercial DNA remover solution to ensure the area is DNA free at the start of your run and re-run the test.
2. If the problem persists then the kit has become contaminated and it will have to be discarded and replaced with a new kit. When you open the new kit, run a simple test to show that changing the kit has solved the problem. Prepare a test which includes only the Positive Control, the Negative Control and one ‘mock sample’. For the ‘mock sample’ add water instead of any sample RNA. The result for the Negative Control and the mock sample should be negative indicating that contamination is no longer present.

### **Preventive action**

An ideal lab set-up has a ‘Clean area’ where the test reagents are prepared and a ‘sample area’ where DNA/RNA samples and the Positive Control template are handled. The best workflow involves setting up all the test components (excluding the positive control template) in the clean area and then moving the tests to the sample area for sample and Positive Control addition. If this method is followed then the kit components are always kept away from possible sources of contamination. For extra security the Negative Control can be completely prepared and sealed in the clean area. All work areas should be decontaminated regularly with DNA remover.

## “Sample preparation failed”

### Explanation

The test has failed because the quality of the sample was not high enough. The Internal Extraction Control component identifies whether the sample has been prepared correctly and is of suitable quality. This error message means that this quality control test has failed and the sample quality is not high enough for analysis.

### Solutions

1. Check the sample preparation protocol for any user errors then repeat.
2. Poor quality samples can result from overloading the sample preparation protocol with too much starting material. Try reducing the amount of starting material then repeat.
3. Failing to add the Internal Extraction Control RNA to your sample during the sample preparation protocol can also lead to a reported result of “sample preparation failed”. Ensure that this step has not been overlooked or forgotten. If your samples are derived from an archive store or from a process separate from your genesig Easy extraction kit; you must add 5µl of Internal Extraction Control RNA into each 0.5ml of your sample to make it suitable for use on the q16.

## “Positive result, poor quality sample”

### Explanation

The test is positive so if you are only interested in obtaining a ‘present or absent’ answer for your sample then your result is reliable. However, the test contains an Internal Extraction Control component that identifies if the sample is of high quality. This quality control test has failed and the sample is not therefore of high enough quality to accurately calculate the exact copy number of RNA present. If you require quantitative information for your sample then proceed with the solutions below.

### Solutions

1. For appropriate solutions, read the “Sample preparation failed” section of this handbook.

## “Test failed”

### Explanation

The test has failed because the Positive Control has not worked. The Positive Control is present to show that all aspects of the test are working correctly together. When this control test fails, the test as a whole is invalidated. This finding indicates that a problem has occurred in the reaction set-up part of the experiment and has nothing to do with sample preparation.

### Solutions

1. Check the entire workflow and test set-up to look for any user errors, then repeat the test e.g. have the right colour pipettes and solutions been used with the correct tubes?
2. Ensure the positive and negative controls are inserted into the correct wells of your q16.
3. A component of the test may have ‘gone off’ due to handling errors, incorrect storage or exceeding the shelf life. When you open a new kit, run a simple test to show that changing the kit has solved the problem. Prepare a test which includes only the Positive Control, the Negative Control and one ‘mock sample’. For the ‘mock sample’ add internal control template instead of any sample RNA. If the Positive Control works, the mock sample will now be called as a negative result.

## “Test failed and is contaminated”

### **Explanation**

The Positive Control is indicating test failure, and the Negative Control is indicating test contamination. Please read the “Test Failed” and “Test contamination” sections of this technical support handbook for a further explanation.

### **Solution**

1. For appropriate solutions, read both the “Test failed” and “Test contaminated” sections of this handbook.

# Sandfly fever Sicilian virus

Sandfly fever Sicilian virus (SFSV) is a distinct arthropod-borne phlebovirus transmitted by sandflies, specifically by *Phlebotomus papatasi*. It was discovered in Sicily, Italy where it affected troops of the World War II Allied Army Forces after the Sicily landings in 1943. The virus causes sandfly fever (also known as Pappataci fever, Phlebotomus fever and three-day fever). It is a single stranded negative- strand RNA virus with a linear genome that is between 11400 and 14700 nucleotides long. The genome comprises 3 segments, termed L , M and S which encode the viral polymerase, structural glycoproteins proteins and nucleocapsid proteins respectively. Virions consist of an envelope and a nucleocapsid and is between 80-100nm in diameter. The virions are spherical and the envelope surrounds three nucleocapsids and has surface projections.

The virus is transmitted by the bite of an infected female sand fly. While *Phlebotomus papatasi* is the main vector, *P. perniciosus* and *P. perfiliewsi* can also transmit the disease. The range of the sandflies determines the spread of the disease and the disease is prevalent in tropical countries of the Eastern hemisphere between latitude 20° and 45° N. The sand fly itself can become infected as a result of biting an infected person any time from 48 hours before until 24 hours after the onset of fever.

In human hosts, the virus multiplies and becomes widely disseminated throughout the body. Within two and one-half to five days after exposure, there is suddenly a feeling of lassitude, abdominal distress, and dizziness, followed within one day by a chilly sensation and a rapid rise in temperature during the next day or two. Symptoms include severe frontal headache and postorbital pain, intense muscular and joint pains, and a flushed appearance of the face. Usually after two days the temperature slowly returns to normal; only rarely is there a second episode of fever. Following the febrile period, there is great fatigue and weakness, accompanied by slow pulse and frequently subnormal blood pressure. Sandfly fever viruses (SFVs) comprise three major serotypes; SFSV, Sandfly fever Naples virus (SFNV) and Toscana virus (TOSV). Sandfly fever Cyprus virus (SFCV), a variant of Sicilian serotype, is also responsible for acute febrile disease and TOSV may cause aseptic meningitis/encephalitis in the infected individuals.

# Specificity

The Primerdesign genesig Kit for Sandfly fever Sicilian virus (SFSV) genomes is designed for the in vitro quantification of SFSV genomes. The kit is designed to have a broad detection profile. Specifically, the primers represent 100% homology with over 95% of the NCBI database reference sequences available at the time of design.

The dynamics of genetic variation means that new sequence information may become available after the initial design. Primerdesign periodically reviews the detection profiles of our kits and when required releases new versions.

If you require further information, or have a specific question about the detection profile of this kit then please send an e.mail to [enquiry@primerdesign.co.uk](mailto:enquiry@primerdesign.co.uk) and our bioinformatics team will answer your question.

## Kit storage and stability

This lyophilised kit is stable at room temperature but should be stored at -20°C on arrival. Once the lyophilised components have been resuspended they should not be exposed to temperatures above -20°C for longer than 30 minutes at a time and unnecessary repeated freeze/thawing should be avoided. The kit is stable for six months from the date of resuspension under these circumstances.

Primerdesign does not recommend using the kit after the expiry date stated on the pack.

## Suitable sample material

All kinds of sample material suited for PCR amplification can be used. Please ensure the samples are suitable in terms of purity, concentration, and RNA/DNA integrity.

## Dynamic range of test

Under optimal PCR conditions genesig SFSV detection kits have very high priming efficiencies of >95% and can detect less than 100 copies of target template.

## Notices and disclaimers

This product is developed, designed and sold for research purposes only. It is not intended for human diagnostic or drug purposes or to be administered to humans unless clearly expressed for that purpose by the Food and Drug Administration in the

USA or the appropriate regulatory authorities in the country of use. During the warranty period Primerdesign genesig detection kits allow precise and reproducible data recovery combined with excellent sensitivity. For data obtained by violation to the general GLP guidelines and the manufacturer's recommendations the right to claim under guarantee is expired. PCR is a proprietary technology covered by several US and foreign patents. These patents are owned by Roche Molecular Systems Inc. and have been sub-licensed by PE Corporation in certain fields. Depending on your specific application you may need a license from Roche or PE to practice PCR. Additional information on purchasing licenses to practice the PCR process may be obtained by contacting the Director of Licensing at Roche Molecular Systems, 1145 Atlantic Avenue, Alameda, CA 94501 or Applied Biosystems business group of the Applied Biosystems Corporation, 850 Lincoln Centre Drive, Foster City, CA 94404. In addition, the 5' nuclease assay and other homogeneous amplification methods used in connection with the PCR process may be covered by U. S. Patents 5,210,015 and 5,487,972, owned by Roche Molecular Systems, Inc, and by U.S. Patent 5,538,848, owned by The Perkin-Elmer Corporation.

## Trademarks

Primerdesign™ is a trademark of Primerdesign Ltd.

genesig® is a registered trademark of Primerdesign Ltd.