

Primerdesign™ Ltd

African swine fever virus

major capsid protein p72 gene

genesig® Standard Kit

150 tests

GENESIG

Kits by Primerdesign

For general laboratory and research use only

Introduction to African swine fever virus

African Swine Fever Virus (ASFV) is a widespread disease which infects members of the pig family (Suidae). A number of tick species are believed to be the vector for the disease, as well as being transmitted by raw pork and pig excrement [1]. After firstly being identified in Kenya in 1921, ASFV became endemic in sub-Saharan Africa, with regular outbreaks being reported across Europe, Asia and South America throughout the century [2]. More recently the virus was introduced in Georgia and spread throughout the region, as well as mass outbreaks occurring in China in 2018 [3].

ASFV is the only member of the Asfaridae family. It is a large enveloped double stranded DNA virus of icosahedral morphology with an average diameter of 200nm and isolates contain genomes between 170-190Kbp encoding for up to 167 open reading frames [2]. The morphology of ASFV consist of several concentric domains. An inner core contains the nucleoid coated with a thick protein layered core shell, which is surrounded by an inner lipid envelope , all of which is encompassed by the capsid [2]. ASFV begins its replication cycle in the nucleus of infected cells before moving to the cytoplasm where the majority of the replication takes place [2]. Gene transcription is highly regulated, with distinct classes of mRNA identified to accumulate at early, intermediate and late transcripts of the virus [2]. The disease induces acute haemorrhagic disease within its hosts, causing high fevers and skin haemorrhages, with death often occurring within ten days of clinical symptoms appearing [4].

References:

- 1: The Centre for Food Security and Public Health (2015), African Swine Fever.
- 2: Galindo, I. and Alonso, C., 2017. African swine fever virus: a review. *Viruses*, 9(5), p.103.
- 3: Zhou, X., Li, N., Luo, Y., Liu, Y., Miao, F., Chen, T., Zhang, S., Cao, P., Li, X., Tian, K. and Qiu, H.J., 2018. Emergence of African swine fever in China, 2018. *Transboundary and emerging diseases*, 65(6), pp.1482-1484.
- 4: Gallardo, C., Ademun, A.R., Nieto, R., Nantima, N., Arias, M., Martín, E., Pelayo, V. and Bishop, R.P., 2011. Genotyping of African swine fever virus (ASFV) isolates associated with disease outbreaks in Uganda in 2007. *African Journal of biotechnology*, 10(17), pp.3488-3497.

Specificity

The Primerdesign genesig kit for African Swine Fever Virus (ASFV) is designed for the in vitro quantification of ASFV genomes. The kit is designed to have a broad detection profile. Specifically, the primers represent 100% homology with over 95% of the NCBI database reference sequences available at the time of design.

The dynamics of genetic variation means that new sequence information may become available after the initial design. Primerdesign periodically reviews the detection profiles of our kits and when required releases new versions.

The target sequence (p72) has been shown to be a good genetic marker for ASFV in other real time PCR based studies (Tignon et al., 2011). This kit was designed to detect recent isolates from the Chinese outbreak of ASFV and, through the use of multiple probes is predicted to detect strains that have been isolated throughout Europe and Africa.

If you require further information, or have a specific question about the detection profile of this kit then please send an e.mail to enquiry@primerdesign.co.uk and our bioinformatics team will answer your question.

Kit contents

- **ASFV specific primer/probe mix (150 reactions BROWN)**
FAM labelled
- **ASFV positive control template (for Standard curve RED)**
- **RNase/DNase free water (WHITE)**
for resuspension of primer/probe mixes
- **Template preparation buffer (YELLOW)**
for resuspension of positive control template and standard curve preparation

Reagents and equipment to be supplied by the user

Real-time PCR Instrument

Extraction kit

This kit is recommended for use with genesig Easy DNA/RNA extraction kit. However, it is designed to work well with all processes that yield high quality RNA and DNA with minimal PCR inhibitors.

oasig™ lyophilised or Precision® PLUS 2X qPCR Master Mix

This kit is intended for use with oasig or PrecisionPLUS 2X qPCR Master Mix.

Pipettors and Tips

Vortex and centrifuge

Thin walled 1.5 ml PCR reaction tubes

Kit storage and stability

This kit is stable at room temperature but should be stored at -20°C on arrival. Once the lyophilised components have been resuspended they should not be exposed to temperatures above -20°C for longer than 30 minutes at a time and unnecessary repeated freeze/thawing should be avoided. The kit is stable for six months from the date of resuspension under these circumstances.

If a standard curve dilution series is prepared this can be stored frozen for an extended period. If you see any degradation in this serial dilution a fresh standard curve can be prepared from the positive control.

Primerdesign does not recommend using the kit after the expiry date stated on the pack.

Suitable sample material

All kinds of sample material suited for PCR amplification can be used. Please ensure the samples are suitable in terms of purity, concentration, and DNA integrity. Always run at least one negative control with the samples. To prepare a negative-control, replace the template DNA sample with RNase/DNase free water.

Dynamic range of test

Under optimal PCR conditions genesig ASFV detection kits have very high priming efficiencies of >95% and can detect less than 100 copies of target template.

Notices and disclaimers

This product is developed, designed and sold for research purposes only. It is not intended for human diagnostic or drug purposes or to be administered to humans unless clearly expressed for that purpose by the Food and Drug Administration in the USA or the appropriate regulatory authorities in the country of use. During the warranty period Primerdesign genesig detection kits allow precise and reproducible data recovery combined with excellent sensitivity. For data obtained by violation to the general GLP guidelines and the manufacturer's recommendations the right to claim under guarantee is expired. PCR is a proprietary technology covered by several US and foreign patents. These patents are owned by Roche Molecular Systems Inc. and have been sub-licensed by PE Corporation in certain fields. Depending on your specific application you may need a license from Roche or PE to practice PCR. Additional information on purchasing licenses to practice the PCR process may be obtained by contacting the Director of Licensing at Roche Molecular Systems, 1145 Atlantic Avenue, Alameda, CA 94501 or Applied Biosystems business group of the Applied Biosystems Corporation, 850 Lincoln Centre Drive, Foster City, CA 94404. In addition, the 5' nuclease assay and other homogeneous amplification methods used in connection with the PCR process may be covered by U.S. Patents 5,210,015 and 5,487,972, owned by Roche Molecular Systems, Inc., and by U.S. Patent 5,538,848, owned by The Perkin-Elmer Corporation.

Trademarks

Primerdesign™ is a trademark of Primerdesign Ltd.

genesig® is a registered trademark of Primerdesign Ltd.

The PCR process is covered by US Patents 4,683,195, and 4,683,202 and foreign equivalents owned by Hoffmann-La Roche AG. BI, ABI PRISM®, GeneAmp® and MicroAmp® are registered trademarks of the Applied Biosystems (Applied Biosystems Corporation). BIOMEK® is a registered trademark of Beckman Instruments, Inc.; iCycler™ is a registered trademark of Bio-Rad Laboratories, Rotor-Gene is a trademark of Corbett Research. LightCycler™ is a registered trademark of the Idaho Technology Inc. GeneAmp®, TaqMan® and AmpliTaqGold® are registered trademarks of Roche Molecular Systems, Inc., The purchase of the Primerdesign reagents cannot be construed as an authorization or implicit license to practice PCR under any patents held by Hoffmann-La Roche Inc.

Principles of the test

Real-time PCR

A ASFV specific primer and probe mix is provided and this can be detected through the FAM channel.

The primer and probe mix provided exploits the so-called TaqMan® principle. During PCR amplification, forward and reverse primers hybridize to the ASFV DNA. A fluorogenic probe is included in the same reaction mixture which consists of a DNA probe labeled with a 5`-dye and a 3`-quencher. During PCR amplification, the probe is cleaved and the reporter dye and quencher are separated. The resulting increase in fluorescence can be detected on a range of qPCR platforms.

Positive control

For copy number determination and as a positive control for the PCR set up, the kit contains a positive control template. This can be used to generate a standard curve of ASFV copy number / Cq value. Alternatively the positive control can be used at a single dilution where full quantitative analysis of the samples is not required. Each time the kit is used, at least one positive control reaction must be included in the run. A positive result indicates that the primers and probes for detecting the target ASFV gene worked properly in that particular experimental scenario. If a negative result is obtained the test results are invalid and must be repeated. Care should be taken to ensure that the positive control does not contaminate any other kit component which would lead to false-positive results. This can be achieved by handling this component in a Post PCR environment. Care should also be taken to avoid cross-contamination of other samples when adding the positive control to the run. This can be avoided by sealing all other samples and negative controls before pipetting the positive control into the positive control well.

Negative control

To validate any positive findings a negative control reaction should be included every time the kit is used. For this reaction the RNase/DNase free water should be used instead of template. A negative result indicates that the reagents have not become contaminated while setting up the run.

Resuspension protocol

To minimize the risk of contamination with foreign DNA, we recommend that all pipetting be performed in a PCR clean environment. Ideally this would be a designated PCR lab or PCR cabinet. Filter tips are recommended for all pipetting steps.

1. Pulse-spin each tube in a centrifuge before opening.

This will ensure lyophilised primer and probe mix is in the base of the tube and is not spilt upon opening the tube.

2. Resuspend the kit components in the RNase/DNase free water supplied, according to the table below.

To ensure complete resuspension, vortex each tube thoroughly.

Component - resuspend in water	Volume
Pre-PCR pack	
ASFV primer/probe mix (BROWN)	165 µl

3. Resuspend the positive control template in the template preparation buffer supplied, according to the table below:

To ensure complete resuspension, vortex the tube thoroughly.

Component - resuspend in template preparation buffer	Volume
Post-PCR heat-sealed foil	
ASFV Positive Control Template (RED) *	500 µl

* This component contains high copy number template and is a VERY significant contamination risk. It must be opened and handled in a separate laboratory environment, away from the other components.

qPCR detection protocol

- 1. For each DNA sample prepare a reaction mix according to the table below:**
Include sufficient reactions for positive and negative controls.

Component	Volume
oasig or PrecisionPLUS 2X qPCR Master Mix	10 μ l
ASFV primer/probe mix (BROWN)	1 μ l
RNase/DNase free water (WHITE)	4 μ l
Final Volume	15 μl

- 2. Pipette 15 μ l of this mix into each well according to your qPCR experimental plate set up.**
- 3. Prepare DNA templates for each of your samples.**
- 4. Pipette 5 μ l of DNA template into each well, according to your experimental plate set up.**
For negative control wells use 5 μ l of RNase/DNase free water. The final volume in each well is 20 μ l.
- 5. If a standard curve is included for quantitative analysis, prepare a reaction mix according to the table below:**

Component	Volume
oasig or PrecisionPLUS 2X qPCR Master Mix	10 μ l
ASFV primer/probe mix (BROWN)	1 μ l
RNase/DNase free water (WHITE)	4 μ l
Final Volume	15 μl

6. Preparation of a standard curve dilution series.

- 1) Pipette 90µl of template preparation buffer into 5 tubes and label 2-6
- 2) Pipette 10µl of Positive Control Template (**RED**) into tube 2
- 3) Vortex thoroughly
- 4) Change pipette tip and pipette 10µl from tube 2 into tube 3
- 5) Vortex thoroughly

Repeat steps 4 and 5 to complete the dilution series

Standard Curve	Copy Number
Tube 1 Positive control (RED)	2×10^5 per µl
Tube 2	2×10^4 per µl
Tube 3	2×10^3 per µl
Tube 4	2×10^2 per µl
Tube 5	20 per µl
Tube 6	2 per µl

7. Pipette 5µl of standard template into each well for the standard curve according to your experimental plate set up.
The final volume in each well is 20µl.

qPCR amplification protocol

Amplification conditions using oasis or PrecisionPLUS 2X qPCR Master Mix.

	Step	Time	Temp
Cycling x50	Enzyme activation	2 min	95 °C
	Denaturation	10 s	95 °C
	DATA COLLECTION *	60 s	60 °C


* Fluorogenic data should be collected during this step through the FAM channel

Interpretation of results


Target	Positive control	Negative control	Interpretation
+	+	-	POSITIVE QUANTITATIVE RESULT calculate copy number
-	+	-	NEGATIVE RESULT
+ / -	+	≤ 35	EXPERIMENT FAILED due to test contamination
+ / -	+	> 35	*
+ / -	-	+ / -	EXPERIMENT FAILED

Positive control template (RED) is expected to amplify between Cq 16 and 23. Failure to satisfy this quality control criterion is a strong indication that the experiment has been compromised

*Where the test sample is positive and the negative control is positive with a Cq > 35 , the sample must be reinterpreted based on the relative signal strength of the two results:

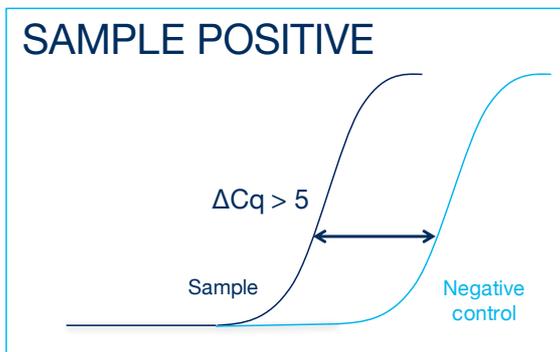

If the sample amplifies > 5 Cq earlier than the negative control then the sample should be reinterpreted (via the table above) with the negative control verified as negative.

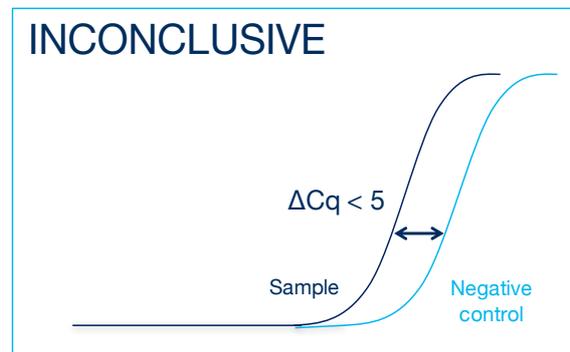

If the sample amplifies < 5 Cq earlier than the negative control then the positive sample result is invalidated and the result should be determined inconclusive due to test contamination. The test for this sample should be repeated.